



610.756.7884


[Home](#) » [Products](#) » [Halsted Mosquito Forceps - Straight, 5"](#)

HALSTED MOSQUITO FORCEPS - STRAIGHT, 5"

Product Number: 17-1450

Halsted Mosquito Forceps, Standard, Straight, 5" Overall Length

Sklar® Halsted Mosquito Forceps are used in multiple procedures as a hemostat for small blood vessels. The ratcheted finger ring handle provides a secure, locking grip to halt blood flow. They are not appropriate for use with large vessels or arteries. A toothed variation is available for grasping delicate tissue in skin grafting, biopsies, or ophthalmologic procedures. Mosquito forceps are known for their fine tips, and short, fully serrated jaws. They are available curved or straight. The standard length for this instrument is 5 inches, although Sklar does offer some with a length of 5-1/2 inches. Halsted Mosquito Forceps are very similar to Hartmann Mosquito Forceps, however, the Hartmann style is finer and more delicate. This product is straight with serrated tips and a length of 5 inches.

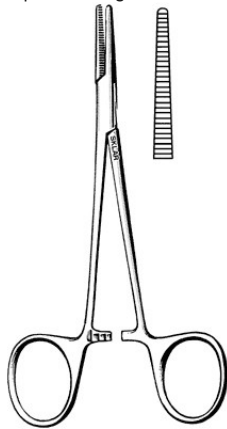

Additional Information

Surname	Halsted Mosquito
Length	5"
Instrument Type	Hemostatic Forceps
Curvature	Straight
Working Surface Style	Serrated
Handle	Finger Ring
Material	Stainless Steel
Disposable or Reusable	Reusable
Sterile or Non-Sterile	Non-Sterile
Latex or Latex-Free	Latex-Free
Grade	Premium OR-Grade

© 2016 Sklar Surgical Instruments. All Rights Reserved.

889 South Matlack Street, West Chester, PA 19382 • Phone: 610.756.7884 • Fax: 610.430.3941 • Service@SklarCorp.com