



610.756.7884


[Home](#) » [Products](#) » [TC Halsey Needle Holder](#)

TC HALSEY NEEDLE HOLDER

Product Number: 21-8001

5" TC Halsey Needle Holder Serrated Tungsten Carbide

Sklar® TC Halsey Needle Holders are ratcheted, finger ring instruments used for holding small needles for suturing, most commonly during plastic surgery and procedures with delicate tissue, such as dermatology. They are available with smooth or delicately serrated jaws. Smooth jaws will cause the least damage to the needle; however, serrated jaws hold the needle more tightly. Physician preference typically determines the type of working surface utilized. Sklar permanently bonds tungsten carbide inserts to the instrument, which maximizes control and precision. This product is straight with cross-serrated tips and a length of 5 inches.

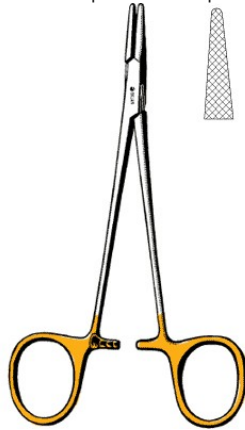

Additional Information

Surname	Halsey / Webster
Length	5"
Tip Configuration	Blunt
Curvature	Straight
Working Surface Style	Cross-Serrated
Handle	Finger Rings / Dual Gold
Material	Tungsten Carbide / Stainless Steel
Disposable or Reusable	Reusable
Sterile or Non-Sterile	Non-Sterile
Latex or Latex-Free	Latex-Free
Grade	Premium OR-Grade

© 2016 Sklar Surgical Instruments. All Rights Reserved.

889 South Matlack Street, West Chester, PA 19382 • Phone: 610.756.7884 • Fax: 610.430.3941 • Service@SklarCorp.com