



610.756.7884


[Home](#) » [Products](#) » [Frazier Suction Tube](#)

## FRAZIER SUCTION TUBE

Product Number: 50-2006

Frazier Suction Tube 6 French, 7-1/2" overall, working length 4" from bend

Sklar's® Frazier Suction Tube is a thin instrument used for removing fluid and debris from confined surgical sites such as the nasal cavity, lumbar and cervical regions or intracranially. Suction is controlled by a small hole on the handle and multiple tip sizes are available. The size of the suction tube is indicated by the French measurement scale. This suction tube also feature a thin wire that can be placed into the tube to remove any tissue that may become trapped throughout the suctioning. This product is a size 6 French, angled tube with a 7-1/2 inches overall, working length 4 inches from bend.

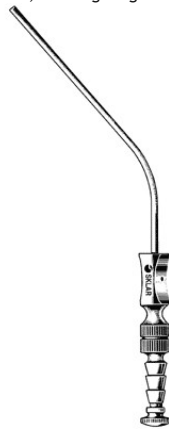

### Additional Information

Surname	Frazier
Length	7-1/2"
Tip Dimensions	6 French
Instrument Type	Suction Tube
Curvature	Curved
Size / Model	working length 4"
Material	Stainless Steel
Disposable or Reusable	Reusable
Sterile or Non-Sterile	Non-Sterile
Latex or Latex-Free	Latex-Free
Grade	Premium OR-Grade

© 2016 Sklar Surgical Instruments. All Rights Reserved.

889 South Matlack Street, West Chester, PA 19382 • Phone: 610.756.7884 • Fax: 610.430.3941 • [Service@SklarCorp.com](mailto:Service@SklarCorp.com)