



610.756.7884


[Home](#) » [Products](#) » [Sklar Edge TC Strabismus Scissors](#)

SKLAR EDGE TC STRABISMUS SCISSORS

Product Number: 16-2705

4-1/2" Sklar Edge Tungsten Carbide Strabismus Scissors Curved Blunt/Blunt

Sklar Edge™ Tungsten Carbide Strabismus Scissors provide an exceptionally durable and hard cutting edge. They are an ophthalmic dissecting instrument and are primarily used during procedures for the treatment of Strabismus, otherwise known as "cross-eyes." These scissors are designed to stay sharp for a significantly longer period of time than standard scissors. Sklar permanently bonds tungsten carbide inserts to the blades of the scissors which maximizes control and precision. Sklar Edge™ scissors can be identified by two gold coated finger rings. This product is curved with smooth, blunt/blunt tips and a length of 4-1/2 inches.

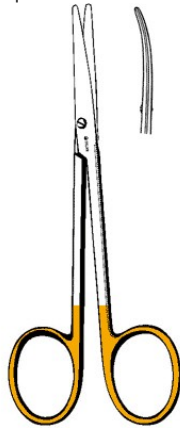

Additional Information

Length	4-1/2"
Tip Configuration	Blunt/Blunt
Instrument Type	Scissors
Curvature	Curved
Working Surface Style	Smooth
Material	Tungsten Carbide / Stainless Steel
Disposable or Reusable	Reusable
Sterile or Non-Sterile	Non-Sterile
Latex or Latex-Free	Latex-Free
Grade	Premium OR-Grade
Trademark	Sklar Edge™ TC

© 2016 Sklar Surgical Instruments. All Rights Reserved.

889 South Matlack Street, West Chester, PA 19382 • Phone: 610.756.7884 • Fax: 610.430.3941 • Service@SklarCorp.com