



610.756.7884


[Home](#) » [Products](#) » [Weitlaner Self-Retaining Retractor](#)

## WEITLANER SELF-RETAINING RETRACTOR

Product Number: 22-8355

5-1/2" Weitlaner Self-Retaining Retractor Blunt 3x4 Teeth

Sklar's® Weitlaner Retractor is a self-retaining, finger ring retractor with a cam ratchet lock used for holding back tissue and exposing a surgical site. It has curved shanks that lead to blades that have between two and six, sharp or blunt, outward curving prongs. The prongs interdigitate when the instrument is closed. The Weitlaner retractor is a popular instrument, most commonly used in basic plastic surgery, large bone and joint procedures, and mastoid surgeries. It is also useful with small, deep incisions and soft tissue dissection at a superficial level. This product has 3x4 blunt prongs and a length of 5-1/2 inches.

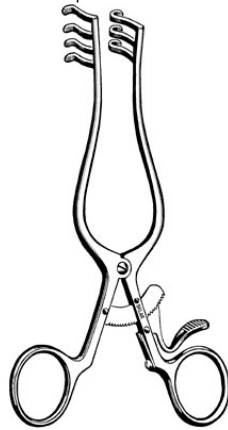

### Additional Information

Surname	Weitlaner
Length	5-1/2"
Tip Configuration	3x4 Teeth
Instrument Type	Retractor
Curvature	Curved
Handle	Finger Rings
Material	Stainless Steel
Disposable or Reusable	Reusable
Sterile or Non-Sterile	Non-Sterile
Latex or Latex-Free	Latex-Free
Grade	Premium OR-Grade

© 2016 Sklar Surgical Instruments. All Rights Reserved.

889 South Matlack Street, West Chester, PA 19382 • Phone: 610.756.7884 • Fax: 610.430.3941 • [Service@SklarCorp.com](mailto:Service@SklarCorp.com)