



610.756.7863


[Home](#) » [Products](#) » [Hartmann Alligator Ear Forceps](#)

HARTMANN ALLIGATOR EAR FORCEPS

Product Number: 67-1045

Hartmann Alligator Ear Forceps Serrated

Sklar's® Hartmann Alligator Ear Forceps are a non-ratcheted, finger ring, ENT instrument commonly used in myringotomy procedures to insert or remove tubes, dressings, or packing. This product has 3-1/2 inch shanks that angle at a screw lock to a 2-1/2 inch shaft and 8mm long triangular tips. The tips are horizontally serrated and each have a central groove to securely grip surgical materials as well as granulation tissue and cerumen.

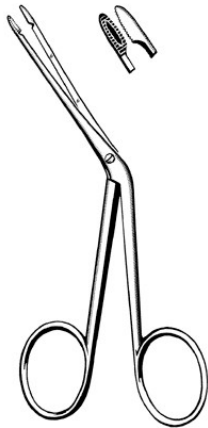

Additional Information

Surname	Hartmann
Instrument Type	Forceps
Material	Stainless Steel
Disposable or Reusable	Reusable
Sterile or Non-Sterile	Non-Sterile
Latex or Latex-Free	Latex-Free
Grade	Premium OR-Grade

© 2016 Sklar Surgical Instruments. All Rights Reserved.

889 South Matlack Street, West Chester, PA 19382 • Phone: 610.756.7863 • Fax: 610.430.3941 • Service@SklarCorp.com