



610.756.7863


[Home](#) » [Products](#) » [Schroeder Tenaculum Forceps - Flat Tips, 10"](#)

## SCHROEDER TENACULUM FORCEPS - FLAT TIPS, 10"

Product Number: 91-1290

10" Schroeder Tenaculum Forceps, Straight, 1x1 Prongs, Braun Style Flat Tips

Sklar® Schroeder-Braun Tenaculum Forceps are a ratcheted, finger ring instrument commonly used in OB/GYN procedures involving manipulation of the cervix or uterus, most commonly D&C, and vaginal or abdominal hysterectomy. They are also used to grasp, hold, and move vessels and arteries during surgical procedures. The tips of the Schroeder-Braun tenaculum forceps have two flattened prong tips that angle inward at 90 degrees. The prongs overlap, providing a strong hold on the tissue, but are also traumatic. This product has straight shanks and an overall length of 10 inches.

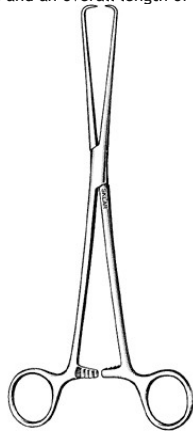

### Additional Information

Surname	Schroeder (sometimes spelled Schroder, Schröder)
Length	10"
Tip Configuration	1x1 Prongs (Braun Style Flat Tips)
Instrument Type	Tenaculum Forceps
Curvature	Straight
Material	Stainless Steel
Disposable or Reusable	Reusable
Sterile or Non-Sterile	Non-Sterile
Latex or Latex-Free	Latex-Free
Grade	Premium OR-Grade

© 2016 Sklar Surgical Instruments. All Rights Reserved.

889 South Matlack Street, West Chester, PA 19382 • Phone: 610.756.7863 • Fax: 610.430.3941 • [Service@SklarCorp.com](mailto:Service@SklarCorp.com)